


CE 0197

**MANI**®

FOR ADVANCED TREATMENT

# PEESO REAMERS GATES DRILLS

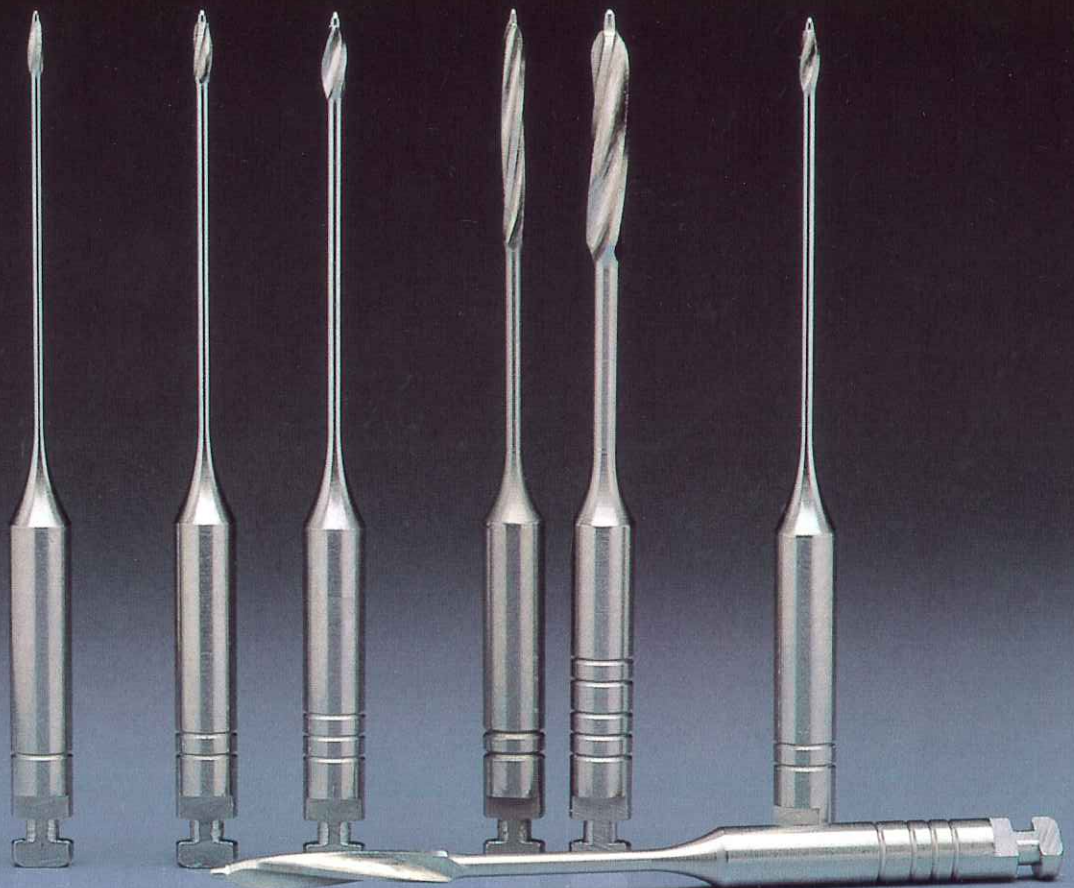

**MANI, INC.**



For advanced treatment ——— We serve with micro processing technology

# PEESO REAMERS GATES DRILLS

**Get  
rusted?**

The material developed by MANI has made it possible to commercialize safe and sharp **PEESO REAMERS & GATES DRILLS** with corrosion resistance to any kind of sterilization including autoclave, chemical vapor, and super acid water.

## The anti-rust material

We adopted "MANI hard-fiber stainless steel" which has uniform and stabilized hardness as the raw material of this instrument.

MANI Hard-Fiber Stainless Steel  
(Austenitic Stainless Steel)




## Length : 28mm, 32mm and 38mm

In addition to standard 32mm length, 28mm and 38mm lengths are available.

You can choose the length exactly you need for molar or canine etc.


### PEESO REAMERS

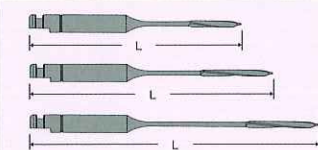

size	Ømm
1	0.70
2	0.90
3	1.10
4	1.30
5	1.50
6	1.70

Length(L) : 28,32,38mm  
Packed : 6pcs. to a package  
Size : 1,2,3,4,5,6 Assorted 1~6

### GATES DRILLS

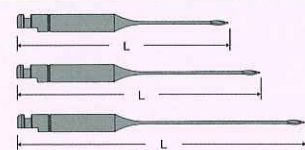

size	Ømm
1	0.50
2	0.70
3	0.90
4	1.10
5	1.30
6	1.50

Length(L) : 28,32,38mm  
Packed : 6pcs. to a package  
Size : 1,2,3,4,5,6 Assorted 1~6

## Safety design

Safety guide of PEESO REAMERS and GATES DRILLS prevent ledges or perforations. Designed to be broken at very near to shank.

## Easy identification

The size of PEESO REAMERS and GATES DRILLS can be easily identified with the number of grooves on the shank.



**Before using the MANI PEESO REAMERS and MANI GATES DRILLS :**

●For dental use only. ●Sterilize PEESO REAMERS and GATES DRILLS before use. ●Insert the shank into the handpiece according to the instruction of the handpiece manufacturer. The insertion should be to full shank length. ●This product should be treated as medical waste when disposed.

**MANI**<sup>®</sup>  
MANI, INC. マニー株式会社

**MANI, INC.**

8-3 KIYOHARA INDUSTRIAL PARK, UTSUNOMIYA, TOCHIGI, 321-3231, JAPAN

Phone : +81-28-667-8592

Fax : +81-28-667-8596

E-mail : dental.exp@ms.mani.co.jp

URL : <http://www.mani.co.jp>

D 00620